



Introducing F&P SleepStyle™

SleepStyle empowers you with the freedom to do what's best for your patients.

We all want the freedom that comes with the ability to choose. To choose how we spend our time. How we go about our work. To choose how we ultimately care for our patients. **With SleepStyle, you can.**





Freedom in simplicity

Get set up and started, simply.

SleepStyle has simplicity woven into its design. Every detail was considered to make it easy for you and your patients to use. In just a few steps, your patients will be on their way to starting therapy.

User-friendly menu.

In collaboration with our users, we created a simple menu. It is designed to save you time in setting up the device and educating patients.



Easy-access chamber.

A simple push of the button gives access to the water chamber which is easy to fill and clean.



Large, responsive buttons.

The layout and function of the buttons work in harmony with the menu to allow quick navigation and ease of use. Patients just need to push "Start".



Built-in connectivity options.

Get information when you and your patients need it. Cellular ensures automated data transfer while **Bluetooth®** technology enables patients to track their progress instantly using the SleepStyle App.



Quiet, integrated design.

The power supply and humidifier have been integrated fully, to minimize space taken on the bedside table.



Powered by technology

Striking the balance between comfort and effective treatment.

The beginning of anything new is usually a crucial time. Learning to adjust to therapy is essential. We've done everything we can to help, from having a world-class auto algorithm with central apnea detection, to providing a full range of innovative comfort options.



ThermoSmart™

Our world-class humidification features a heated breathing tube with AirSpiral™ technology that is designed to minimize condensation. Our enhanced algorithm delivers a comfortable humidification experience for patients.**



SensAwake™

A responsive ramp throughout the night. When wakefulness is detected, treatment pressure is automatically and promptly reduced to a more comfortable level. This eases the return to sleep and helps effective treatment to continue.




Expiratory Relief

An added layer of comfort. Automatically relieves CPAP pressure during exhalation to make breathing more comfortable.

You're in control

SleepStyle is so much more than just a device. Our patient management and support tools complete a seamless experience to help patients succeed in embracing therapy.

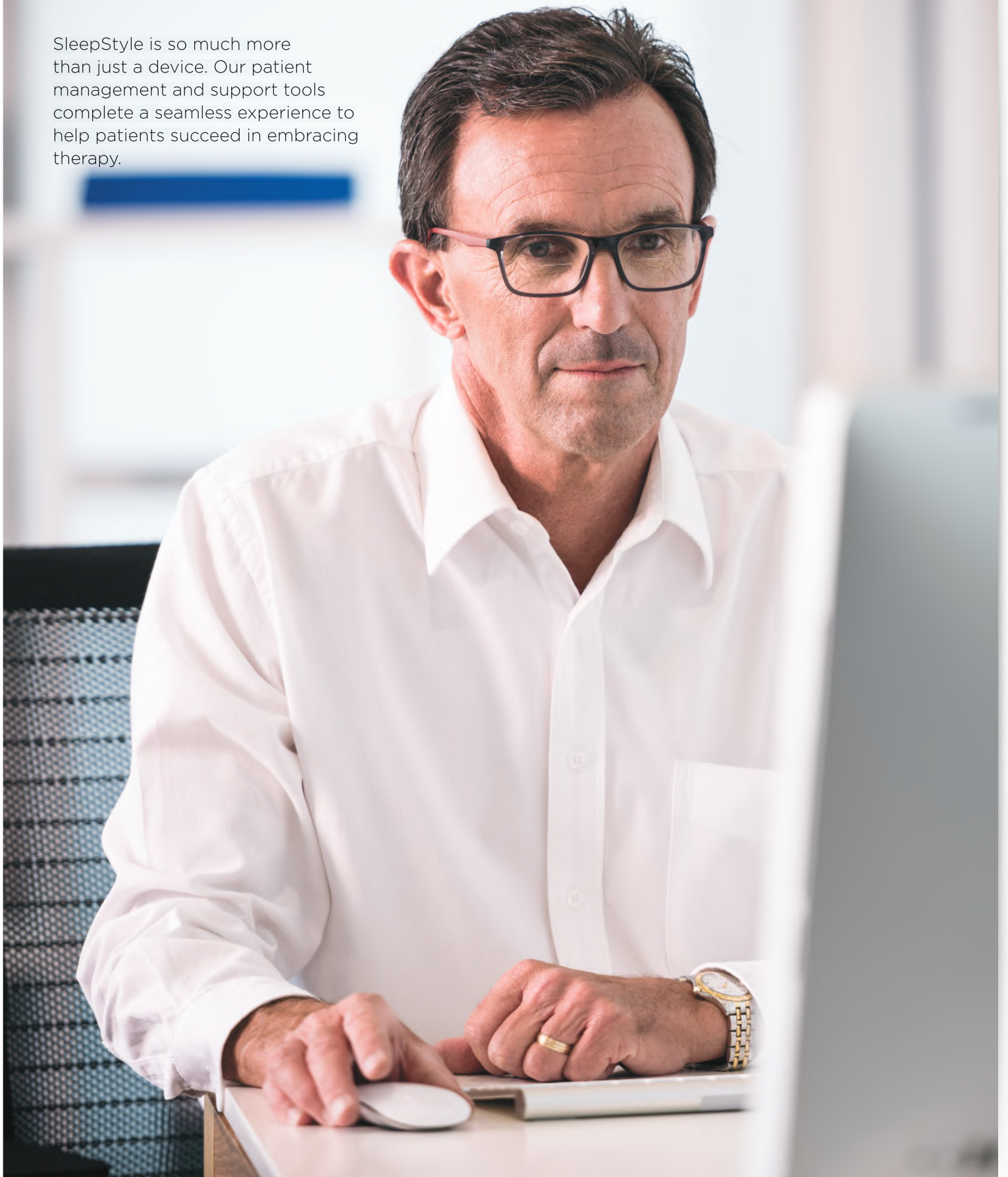

Proactively care for your patients.

F&P InfoSmart™

F&P InfoSmart gives you the freedom to focus your time where it's needed the most. It automatically monitors your patients on an ongoing basis, alerting you to those who need your attention. This will enable you to maximize the chances of every patient's success.



Keep on track with patient involvement.

F&P SleepStyle™ App and Web

Designed to enhance engagement with therapy. Patients have access to meaningful information allowing them to follow their progress and compare it against how they did on previous nights.

Patients can also take action with troubleshooting tips, making the crucial acclimatization period as easy as possible.



Product Specifications

Dimensions	144 H x 177 W x 183 D mm (5.7 H x 7.0 W x 7.2 D in)		
Weight	1.7 kg (3.7 lb) Packed weight (max.): 2.7 kg (5.9 lb)		
Performance	Pressure Range: 4 to 20 cmH ₂ O		
Noise Level	28.0 ± 1.5 dBA		
Water Chamber Volume	380 ml up to the maximum water level line		
Standards Compliance	IEC 60601-1 (Ed 3.1):2012; ISO 80601-2-70:2015; ISO 8185:2007; IEC 60601-1-11:2015; IEC 60601-1-2:2014; Radio Equipment Directive (RED)		
Data Recording	The InfoUSB will store up to 5 years of summary efficacy data, 365 days of detailed efficacy data and 140 hours of high-resolution pressure, leak and flow data. Without an InfoUSB, the device's internal memory is capable of storing up to 1 year of summary efficacy data, 30 days of detailed efficacy data and 20 hours of high-resolution pressure, leak and flow data.		
Electrical Ratings	Rated supply voltage	Rated current input	Rated supply frequency
	100-115 V	1.2 A (2.5 A max.)	50-60 Hz
	220-240 V	1.10 A (2.3 A max.)	50-60 Hz

F&P SleepStyle AUTO/CPAP Performance Features

Performance Features	SleepStyle Auto	SleepStyle CPAP
Fully integrated humidifier	✓	✓
Auto-adjusting pressure	✓	No
ThermoSmart technology**	✓	✓
SensAwake	✓	✓
Expiratory relief	✓	✓
Central Sleep Apnea detection	✓	✓
Ramp	✓	✓
Auto-altitude adjustment	✓	✓
Leak compensation	✓	✓
Efficacy reporting	✓	✓
Compliance reporting	✓	✓
F&P InfoUSB	✓	✓
Bluetooth* wireless technology	✓	✓
Cellular modem*	✓	✓

Compatible with	SleepStyle Auto	SleepStyle CPAP
F&P InfoSmart™	✓	✓
F&P SleepStyle App and Web†	✓	✓

Device Model Product Codes REF

Region	SleepStyle Auto	SleepStyle CPAP
UK/Ireland	SPSAAK*/SPSABK##	SPSCAK*/SPSCBK##
Italy	SPSAAW*/SPSABW##	SPSCAW*/SPSCBW##
Belgium, Netherlands, Luxembourg	SPSAAW*/SPSABW##	SPSCAW*/SPSCBW##

*Not available in all models.

**Not available in all models. A ThermoSmart Heated Breathing Tube is required to activate ThermoSmart.

†Not Available in all countries.

‡Applied Parts – to fit 22 mm (0.86 in.) Conical Connector.

*With modem

##Without modem

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG Inc. and any use of such marks by Fisher & Paykel Healthcare is under license. Other trademarks and trade names are those of their respective owners.

SleepStyle, InfoSmart, ThermoSmart, AirSpiral and SensAwake are trademarks of Fisher & Paykel Healthcare Limited. For patent information, see www.fphcare.com/ip

CE 0123 This medical device is a regulated health product which, according to this regulation, bears the CE mark.

614134 A © 2018 Fisher & Paykel Healthcare Limited

Included with the F&P SleepStyle AUTO/CPAP

1 x Carry Bag
 1 x SleepStyle Device
 1 x Breathing Tube
 1 x Power Cord
 1 x Use and Care Guide
 1 x Water Chamber
 1 x Chamber Seal
 1 x Outlet Seal
 1 x F&P InfoUSB (already in InfoUSB port)
 1 x Air Filter (already in device)
 1 x Spare Air Filter
 1 x Spare Elbow (for use with a standard breathing tube)
 1 x SleepStyle Quick Reference Guide

Replacement Parts

900SPS100 SleepStyle Water Chamber
 900SPS101 SleepStyle Chamber Seal
 900SPS110 SleepStyle Air Filters (2-pack)
 900SPS111 SleepStyle Air Filter Single
 900SPS120 ThermoSmart Breathing Tube††
 900SPS121 Standard Breathing Tube with Elbow††
 900SPS122 SleepStyle Elbow (for use with a Standard Breathing Tube)
 900SPS140 SleepStyle Device Lid
 900SPS141 SleepStyle Outlet Seal
 900SPS142 SleepStyle Carry Bag
 900SPS162 SleepStyle Power Cord – Europe
 900SPS165 SleepStyle Power Cord – UK
 900SW101 F&P InfoUSB

Replacement User Guides

900SPS171 Use and Care Guide – ROW
 900SPS176 Use and Care Guide – East Europe
 900SPS177 Use and Care Guide – West Europe
 900SPS173 SleepStyle Quick Reference Guide